
FRONT

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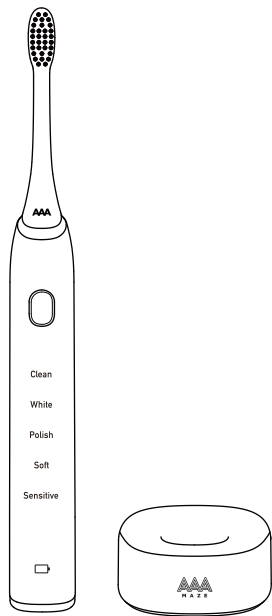ISTRUZIONI D'USO

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
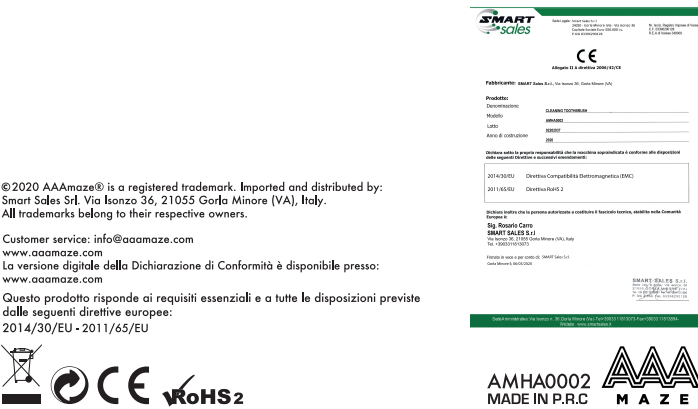
www.aamaze.com - info@aamaze.com

BACK

DYNAMIC SONIC CLEANING TOOTHBRUSH

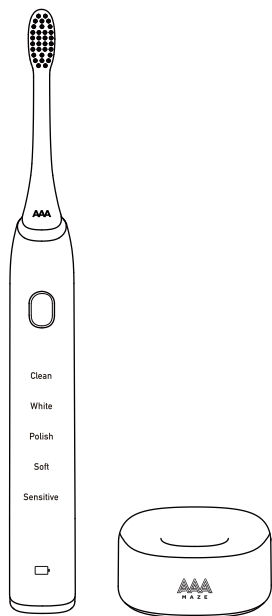INSTRUCTIONS

WHAT INSIDE THE PACK:

1. 1x Dynamic Sonic Cleaning Toothbrush
2. 2x Brush Head
1. 1x Wireless Charging Base
4. Quick Guide

TECHNICAL SPECIFICATIONS:

- PRODUCT SIZE: -Main device : 27.8x27.8x248.6mm
-Charging base: 65.6x65.6x26.3mm
- PRODUCT WEIGHT: -Main body: 115g
-Charging Base: 100g

- COLOR: White
- CHARGING TYPE: Wireless Charging base
- CHARGING TIME: 8 Hours
- BATTERY: 3.7V-2.0W Lithium rechargeable Battery
- CAPACITY: 800mAh
- BATTERY LIFE: 15-20 Days
- BRUSH MODES: 5 different modes, each mode will be lighted up when relevant mode function is on

Clean mode: 34800 strokes per minute
Whiten mode: 33000-40000 strokes per minute
Soft mode: 33000 strokes per minute
Polish Mode: 35000 strokes per minute
Sensitive mode: 28200 strokes per minute

DANGER:

Keep the charging base away from water. Do not place or store it over near water contained in a bathtub, washbasin, etc. Do not immerse the charging base in water or any other liquid. After cleaning make sure the charging base is completely dry before connect it to the mains.

WARNING:

Read all instructions carefully before use and save it for future reference.

- Check the USB cord before you connect the local mains.
- Never use the charging base if it is damaged in any way. If the USB cable is damaged, always have it replaced with one of the original type in order to avoid a hazard.

EN-1

- If the product is damaged in any way (included brush head, toothbrush handle), stop using it immediately.
- Do not use the charging base near heated surface and keep the product away from strong electromagnetic field when charging.
- This product can only be used by children aged from 18 years old and above and by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instructions concerning use of the appliance in a safe way and if they understand the hazard involved. Cleaning and user maintenance shall not be made by children unless they are older than 10 and supervised. Keep the appliance and its cable out of reach of children aged less than 10 years.
- Children shall not play with the product.
- If occur tickling or slight bleeding upon first use, this is normal status. It's easy after first days and disappear absolutely within 1 or 2 weeks.

CAUTION:

- Do not clean the brush head, handle and charging base in the dishwasher. Do not use the toothbrush in the bath or shower.

- If you had oral or gum surgery in the previous 2 months, consult your dentist before using it.

- If you have medical concerns, consult your doctor before you use the appliance. The appliance complies with the safety standard for electromagnetic devices. If you have a pacemaker or other implanted device, consult your doctor before use.

- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose.

- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of bristles fade appear.

- Do not use other brush heads instead of the ones recommended by the manufacturer.

- Supervise children while they are brushing until they have established good brushing habits.

EN-2

MEDICAL WARNINGS:


- Consult your dental professional before you use this product if you have had oral or gum surgery in the previous two months.

- Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after one week of use.

- This product complies with the safety standards for electromagnetic devices.

If you have a pacemaker or other implanted device, contact your physician the device manufacturer prior to use.

PRODUCT COMPONENTS:

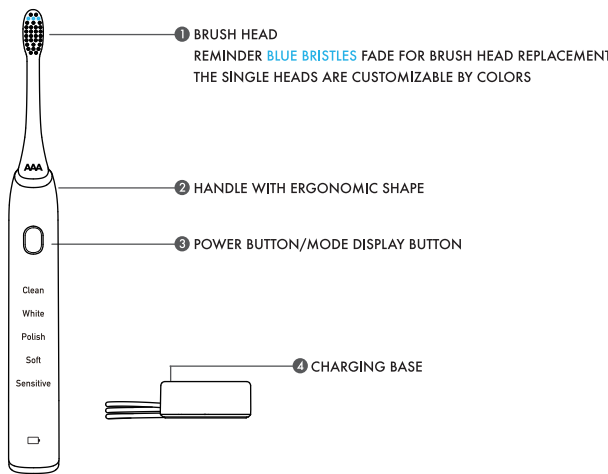EN-3

BATTERY CHARGING:

- Charge for above 12hrs before first use.

- When the battery power is low, the device will vibrate for 4 times and the LED indicator will flash in red to remind the user to charge the battery.

- Put the device into charging base and check the correct connection to the power source. The LED indicator flash in white during battery charging.

- The LED indicator light in white color when battery is fully charged.

BRUSH MODES:

- **CLEAN**
The best for the plaque removal
- **WHITE**
Use the 2-minute whitening mode to remove surface stains, plus an additional 30 seconds to brighten and polish your front teeth
- **POLISH**
Dental polishing, whiten the teeth
- **SOFT**
Extra gentle mode for sensitive teeth and gums
- **SENSITIVE**
Gentle teeth and gum cleaning

34.800
33.000-40.000
35.000
33.000
28.200


EN-4

HOW TO USE IT:

- 01 Put brush head into toothbrush body and leave a small space between head and body
- 02 Take some toothpaste on the wet bristle
- 03 Smudge the toothpaste to teeth before use (Note: launch toothbrush after put into mouth)
- 04 Change the cleaning area based on briefly pause during the brushing cycle (Note: no need to left and right or upper and down to move after launched)

HOW TO CLEAN IT:

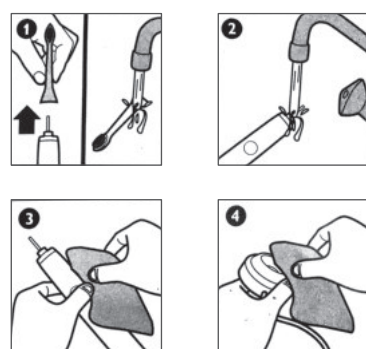Clean the smart toothbrush periodically. Rinse the handle and the brush head. Wipe down all surfaces with a soft cloth and mild soap, if necessary. Do not use sharp objects to clean. Do not clean the brush head, handle, charging base in the dishwasher.

EN-5